

# Advanced Wound Closure

Innovative new ways to close surgical incisions



**Smarter  
wound care  
starts here**

# Designed with doctors to meet today's needs


Now more than ever, you need efficient, effective options that answer all of today's needs in one neat solution. We offer a comprehensive selection of **innovative wound-closure alternatives** designed to **increase efficiency, promote healing, and elevate traditional wound care.**

Ideal for use in the office and the OR, these topical applications provide a smarter, faster way to work while keeping patient comfort and convenience in mind.

A complete  
portfolio of  
wound-closure  
solutions


## Long wounds and incisions such as:

- Shoulders
- C-sections
- Hips
- Knees

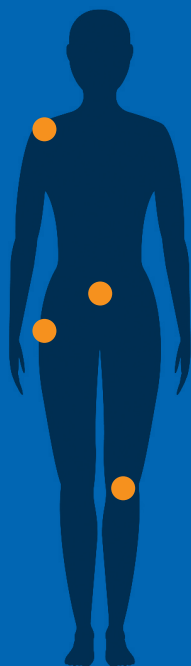

## Short wounds and incisions such as:

- Thyroid
- Pacemaker
- Catheter
- Laparoscopic
- Hernia repair
- Excisions
- Lacerations

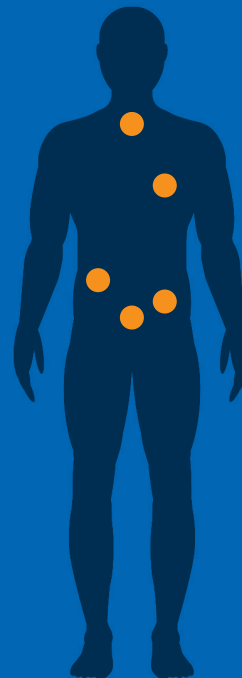


# BandGrip<sup>®</sup>

## Micro-anchor skin closures

Closing long incisions requires a more specialized approach than traditional topical skin adhesive—especially in **orthopedic areas like knees, shoulders and hips**. That's why BandGrip is unlike any other product on the market.

Made with unique, **curved micro-anchors to gently grip the skin**, BandGrip holds surgical incisions and lacerations closed during healing. Designed to avoid common problems associated with conventional closure methods, it can **improve the experience** for patients and caregivers.

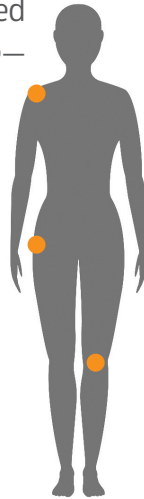


- **Simple, intuitive application** can be successfully implemented by a wide range of caregivers—including the patient
- **Transparent material** makes wound monitoring easier
- **Perforations** allow for wound drainage
- **Minimally invasive design** closes skin without needles or staples, reducing risk of needle sticks and negative cosmetic outcomes
- **Smooth, shower-safe surface** won't snag on dressings or clothes
- **Not made with natural rubber latex**

Made for  
tele-  
medicine?

BandGrip's transparent material makes it ideal for post-op telemedicine. Not only can surgeons monitor patient healing remotely such as on video call, they can easily coach patients through BandGrip changes and removal.

- **Easy application** requires no patient training
- **No need for patients** to return to the facility
- **Frees up surgeon time** for more surgeries vs. post-op patient visits

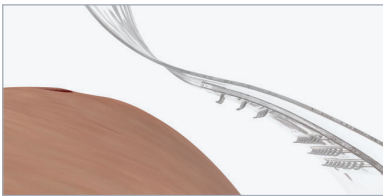


Micro-anchors



# LiquiBand® Rapid™

## Topical skin adhesive


Working smarter and faster in the OR is critical—and easier than ever with LiquiBand Rapid. While other products dry in 90 seconds or more,<sup>1</sup> **LiquiBand Rapid's dry time is as fast as 30 seconds**—which can help **improve OR turnover time** and help you care for more patients, more effectively.

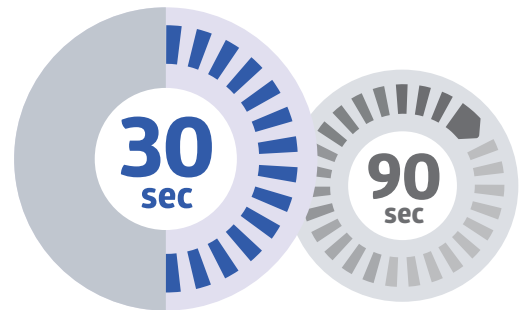


LiquiBand Rapid

Other products

- **More than just fast**—it provides a strong, effective microbial barrier—in vitro studies show no detectable micro-organism penetration<sup>2</sup>
- **Porous foam tip** helps control the application of the adhesive without dripping
- **Accelerant-coated tip** speeds up the adhesive drying process and helps reduce patient turnover time
- **Winged applicator** protects fingers from the risk of glass shards and makes one-handed activation safe and easy to control<sup>3</sup>
- **Highlighter shape** allows for accurate application of broad or thin lines

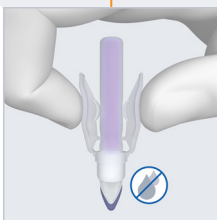

Winged applicator

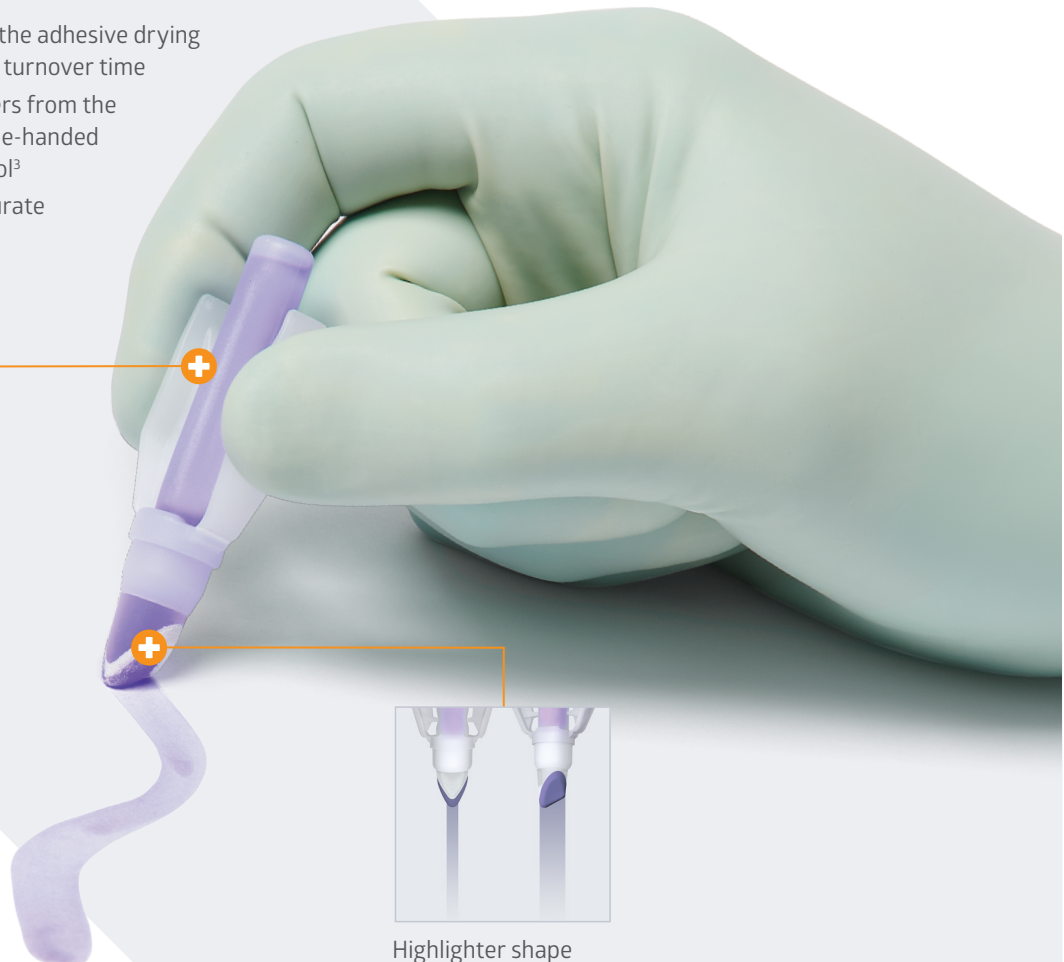


Highlighter shape

# Skin Affix™


## Topical skin adhesive

Skin Affix is an easy-to-use skin adhesive created to close wounds from surgical incisions, **minimally invasive surgical procedures** and simple, thoroughly cleansed trauma-induced lacerations.

- **Low exothermic reaction** means patients do not feel a burning sensation
- **Extended application time** for use on multiple incisions
- **Closes minimally** invasive surgical wounds more quickly than sutures—secure wound closure within 2 minutes
- **High viscosity** allows for more control during application
- **Versatile foam tip** allows for broad and narrow strokes
- **Transparent applicator** for controlled delivery

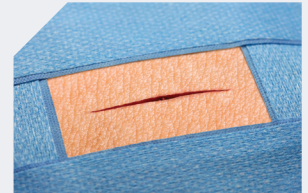

Before

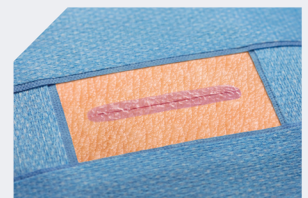

After

### The benefits of 2-octyl cyanoacrylate



2-Octyl cyanoacrylate formulation makes a stronger, flexible film that provides a more durable, long-lasting wound closure.<sup>4</sup>



Can be used in place of sutures that are 5-0 or smaller in diameter.



No need for removal; sloughs off naturally within 5 to 10 days.



# How to order



## BandGrip® Micro-Anchor Skin Closures

Item No.	Description	Pkg.
MSCBG1002	No Perforation, 4 Rows of Anchors, 9 cm x 4 cm	36/bx
MSCBG1002F	With Perforation, 4 Rows of Anchors, 9 cm x 4 cm	36/bx
MSCBG1003	No Perforation, 2 Rows of Anchors, 6.5 cm x 2 cm	36/bx
MSCBG1003F	With Perforation, 2 Rows of Anchors, 6.5 cm x 2 cm	36/bx
MSCBG1005	No Perforation, 2 Rows of Anchors, 5 cm x 7.5 cm	36/bx
MSCBG1005F	With Perforation, 2 Rows of Anchors, 5 cm x 7.5 cm	36/bx



## LiquiBand® Rapid™ Topical Skin Adhesive

Item No.	Description	Pkg.
MSC091050	0.8 ml	36/cs
MSC091050H	0.8 ml	1 ea
MSC091050Z	0.8 ml	6/bx



## Skin Affix™ Topical Skin Adhesive

Item No.	Description	Pkg.
MSC091040	0.55 ml	72/cs
MSC091040H	0.55 ml	1 ea
MSC091040Z	0.55 ml	12/bx

Get the wound-closure coverage your facility needs. To learn more, contact your Medline Representative or visit [medline.com/fasterwoundclosure/](https://medline.com/fasterwoundclosure/)



**We make  
healthcare  
run better™**

**Medline Industries, Inc.**  
Three Lakes Drive, Northfield, IL 60093  
Medline United States | 1-800-MEDLINE (633-5463)  
[medline.com](https://medline.com) | [info@medline.com](mailto:info@medline.com)

**Medline Canada**  
1-800-268-2848 | [medline.ca](https://medline.ca) | [canada@medline.com](mailto:canada@medline.com)  
**Medline México**  
01-800-831-0898 | [medlinemexico.com](https://medlinemexico.com) | [mexico@medline.com](mailto:mexico@medline.com)

Follow us    